


# mRNA-based COVID-19 vaccine boosters induce neutralizing immunity against SARS-CoV-2 Omicron variant

## Graphical abstract

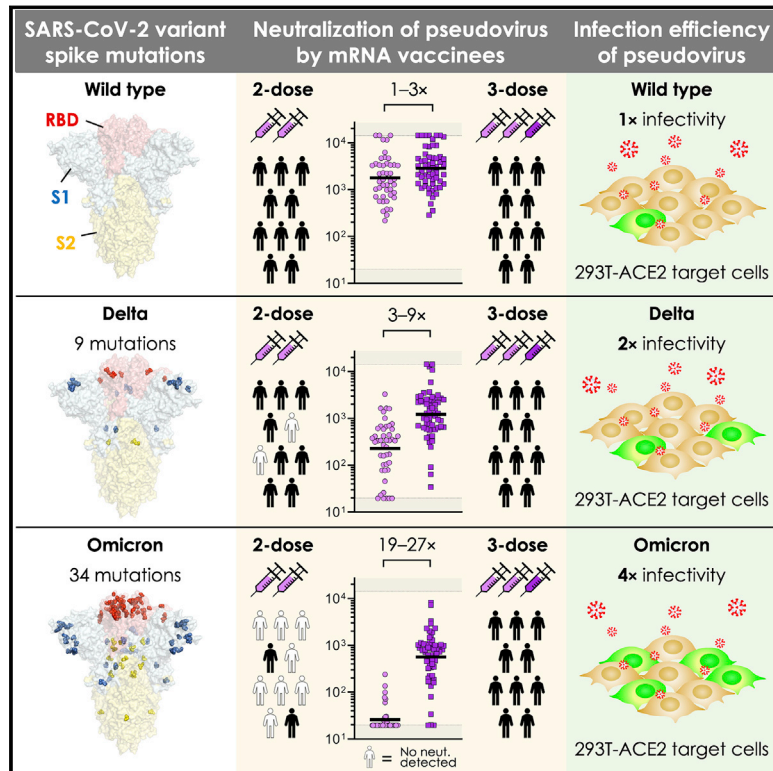

## Authors

Wilfredo F. Garcia-Beltran,  
 Kerri J. St. Denis,  
 Angelique Hoelzemer, ..., A. John Iafrate,  
 Vivek Naranbhai, Alejandro B. Balazs

## Correspondence

wgarcibeltran@mgh.harvard.edu  
 (W.F.G.-B.),  
 abalazs@mgh.harvard.edu (A.B.B.)

## In brief

SARS-CoV-2 Omicron variant pseudovirus exhibits escape from vaccine-induced humoral immunity. However, a third dose of COVID-19 mRNA vaccine elicited humoral immunity capable of cross-neutralizing this strain. In addition, pseudovirus produced with the Omicron spike exhibited more efficient transduction of ACE2-expressing target cells than other variants.

## Highlights

- The SARS-CoV-2 Omicron variant harbors 34 mutations in the spike, more than other variants
- Two doses of mRNA-based vaccines elicit poor neutralization of Omicron
- Three mRNA vaccine doses elicit potent variant cross-neutralization, including Omicron
- The Omicron pseudovirus infects cells more efficiently than other SARS-CoV-2 variants

Article

# mRNA-based COVID-19 vaccine boosters induce neutralizing immunity against SARS-CoV-2 Omicron variant

Wilfredo F. Garcia-Beltran,<sup>1,2,13,\*</sup> Kerri J. St. Denis,<sup>2,13</sup> Angelique Hoelzemer,<sup>3,4,5</sup> Evan C. Lam,<sup>2</sup> Adam D. Nitido,<sup>2,6</sup> Maegan L. Sheehan,<sup>2</sup> Cristhian Berrios,<sup>1</sup> Onosereme Ofoman,<sup>1</sup> Christina C. Chang,<sup>7,8</sup> Blake M. Hauser,<sup>2</sup> Jared Feldman,<sup>2,6</sup> Alex L. Roederer,<sup>2,6</sup> David J. Gregory,<sup>9,10</sup> Mark C. Poznansky,<sup>9,11</sup> Aaron G. Schmidt,<sup>2</sup> A. John Iafrate,<sup>1</sup> Vivek Naranbhai,<sup>8,11,12</sup> and Alejandro B. Balazs<sup>2,14,\*</sup>

<sup>1</sup>Department of Pathology, Massachusetts General Hospital, Boston, MA 02114, USA

<sup>2</sup>Ragon Institute of MGH, MIT, and Harvard, Cambridge, MA 02139, USA

<sup>3</sup>First Department of Internal Medicine, Division of Infectious Diseases, University Medical Centre Eppendorf, Hamburg, Germany

<sup>4</sup>German Center for Infection Research (DZIF), Site Hamburg-Lübeck-Borstel-Riems, Germany

<sup>5</sup>Research Department Virus Immunology, Leibniz Institute for Experimental Virology, Hamburg, Germany

<sup>6</sup>Ph.D. Program in Virology, Division of Medical Sciences, Harvard Medical School, Boston, MA 02115, USA

<sup>7</sup>Alfred Hospital, Central Clinical School, Monash University, Victoria 3181, Australia

<sup>8</sup>Center for the AIDS Programme of Research in South Africa, Durban 4001, South Africa

<sup>9</sup>Vaccine and Immunotherapy Center, Massachusetts General Hospital, Boston, MA 02129, USA

<sup>10</sup>Pediatric Infectious Disease, Massachusetts General Hospital for Children, Boston, MA 02114, USA

<sup>11</sup>Massachusetts General Hospital Cancer Center, Boston, MA 02114, USA

<sup>12</sup>Dana-Farber Cancer Institute, Boston, MA 02215, USA

<sup>13</sup>These authors contributed equally

<sup>14</sup>Lead contact

\*Correspondence: [wgarcibeltran@mgh.harvard.edu](mailto:wgarcibeltran@mgh.harvard.edu) (W.F.G.-B.), [abalazs@mgh.harvard.edu](mailto:abalazs@mgh.harvard.edu) (A.B.B.)  
<https://doi.org/10.1016/j.cell.2021.12.033>

## SUMMARY

Recent surveillance has revealed the emergence of the SARS-CoV-2 Omicron variant (BA.1/B.1.1.529) harboring up to 36 mutations in spike protein, the target of neutralizing antibodies. Given its potential to escape vaccine-induced humoral immunity, we measured the neutralization potency of sera from 88 mRNA-1273, 111 BNT162b, and 40 Ad26.COVS vaccine recipients against wild-type, Delta, and Omicron SARS-CoV-2 pseudoviruses. We included individuals that received their primary series recently (<3 months), distantly (6–12 months), or an additional “booster” dose, while accounting for prior SARS-CoV-2 infection. Remarkably, neutralization of Omicron was undetectable in most vaccinees. However, individuals boosted with mRNA vaccines exhibited potent neutralization of Omicron, only 4–6-fold lower than wild type, suggesting enhanced cross-reactivity of neutralizing antibody responses. In addition, we find that Omicron pseudovirus infects more efficiently than other variants tested. Overall, this study highlights the importance of additional mRNA doses to broaden neutralizing antibody responses against highly divergent SARS-CoV-2 variants.

## INTRODUCTION

The SARS-CoV-2 Omicron variant (BA.1/B.1.1.529) was first detected in Botswana and reported to the World Health Organization (WHO) in November 2021 as a novel variant with an unprecedented number of previously described and novel mutations with immunevasive potential. A subsequent and rapid increase in Omicron cases in South Africa resulted in its designation as a novel variant of concern (VOC) by the WHO ([www.who.int](http://www.who.int)). This variant harbors up to 59 mutations throughout its genome, with as many as 36 of these occurring within the spike protein, the mediator of host cell entry and the main target of neutralizing antibodies. Studies of previous SARS-CoV-2 variants have

demonstrated that mutations within the receptor binding domain (RBD) mediate escape from vaccine-induced neutralizing antibodies (Cele et al., 2021a; Garcia-Beltran et al., 2021a; Zhou et al., 2021), and in some cases, increase infectivity through enhanced affinity for ACE2 (Tian et al., 2021). Omicron RBD contains 15 mutations, some of which overlap with previously studied variants. For example, Beta (B.1.351) and Gamma (P.1) harbor mutations in residues K417, E484, and N501 that potentially diminish vaccine-induced neutralization (Garcia-Beltran et al., 2021a), possibly the result of neutralizing antibody responses being focused toward a limited set of RBD epitopes, as has been previously described (Barnes et al., 2020; Greaney et al., 2021a).

In the United States, three vaccines have been approved by the FDA or are under emergency use authorization (EUA), all of which use the original wild-type SARS-CoV-2 spike protein first identified in Wuhan, China, as the sole immunogen. These are formulated as spike-encoding mRNA in lipid nanoparticles (BNT162b2 manufactured by Pfizer-BioNTech and mRNA-1273 manufactured by Moderna) or as an adenovirus vectored vaccine (Ad26.COV2.S manufactured by Janssen/Johnson & Johnson) (Baden et al., 2021; Polack et al., 2020; Sadoff et al., 2021). These SARS-CoV-2 vaccines have been remarkably successful in inducing neutralizing humoral and cellular immunity and, more importantly, reducing COVID-19 infections, hospitalizations, and deaths in clinical trials (Baden et al., 2021; Polack et al., 2020; Sadoff et al., 2021) and during rapid worldwide deployment (Tregoning et al., 2021). However, it has now been shown that neutralizing antibody responses and vaccine effectiveness vary by vaccine agent, decrease with increased time post vaccination, and are negatively impacted by emerging variants (Bajema et al., 2021; Bar-On et al., 2021; Cromer et al., 2021; Khoury et al., 2021; Lopez Bernal et al., 2021; Naranbhai et al., 2021a; Tregoning et al., 2021). In an effort to combat waning antibody responses and the emergence of new variants, a third dose of mRNA vaccine (“boosters”) has been approved for individuals vaccinated >6 months ago and has been shown to be very effective at inducing high neutralizing antibody titers (Bar-On et al., 2021). In the case of Ad26.COV2.S vaccinees, cross-over mRNA doses were recommended >2 months from primary vaccination. However, while neutralization of wild-type SARS-CoV-2 has been shown to predict the effectiveness of vaccines against variants (Earle et al., 2021; Khoury et al., 2021), it is unclear whether this correlation will be maintained among boosted individuals and for highly mutated variants like Omicron (Bajema et al., 2021; Lopez Bernal et al., 2021; Tregoning et al., 2021).

We previously developed and validated a high-throughput pseudovirus neutralization assay to understand differences in immunity by vaccine and host characteristics against SARS-CoV-2 variants and other coronaviruses (Garcia-Beltran et al., 2021a, 2021b; Naranbhai et al., 2021a, 2021b). Here, we used this assay to test sera from 239 individuals who had been fully vaccinated with one of three vaccines approved in the United States—mRNA-1273, BNT162b2, or Ad26.COV2.S—against wild-type, Delta, and Omicron SARS-CoV-2 pseudoviruses. This included 70 individuals who had received a third dose of mRNA vaccine after primary vaccination series as either a cross-over or solely mRNA-based vaccination regimen. Remarkably, we found that all three primary vaccine series resulted in low-to-absent neutralization of SARS-CoV-2 Omicron. However, mRNA vaccine recipients who received a third dose exhibited potent neutralization against Omicron, despite exhibiting wild-type neutralization titers similar to those in two-dose vaccinated (“non-boosted”) individuals. In addition, *in vitro* infection experiments demonstrated that the Omicron pseudovirus continues to rely upon the human ACE2 receptor for target cell entry and infects target cells 4-fold more efficiently than wild-type pseudovirus and 2-fold more efficiently than Delta pseudovirus. Together, our results highlight that the SARS-CoV-2 Omicron variant evades vaccine-induced neutralizing immunity under current vaccine regimens and is

more infectious than previous variants. Nonetheless, our finding of potent cross-neutralizing immunity against Omicron in individuals that received a third dose of mRNA vaccine suggests that existing vaccines may overcome evasion of humoral immunity by future variants of concern.

## RESULTS


### Emergence of SARS-CoV-2 Omicron as a novel and highly mutated VOC

Over the course of more than 270 million confirmed SARS-CoV-2 infections worldwide, the virus has undergone remarkable diversification, producing >1,500 uniquely identified Pango lineages (Rambaut et al., 2020) ([coronavirus.jhu.edu](https://coronavirus.jhu.edu)). Some of these have demonstrated evidence of increased transmissibility, virulence, and/or immune evasion, prompting the WHO to classify five lineages as current VOCs ([www.who.int](http://www.who.int)). The Omicron variant, also known as PANGO lineage BA.1 or B.1.1.529, was first reported in November 2021 and received its VOC designation within days on account of its unique mutational profile and the dramatic rise in cases observed in Gauteng, South Africa. While the Delta variant is now the dominant SARS-CoV-2 variant worldwide after overtaking the Alpha variant in July 2021, the rise of Omicron infections in regions where Delta is circulating suggests that Omicron may overtake Delta to become the next dominant strain. Despite the substantial recent expansion of the Delta lineage, phylogenetic analysis suggests that the Omicron variant was derived from the Alpha lineage and only recently detected by genomic surveillance (Figure 1A). In comparison to the nine mutations or deletions found in Delta, the Omicron lineage we tested harbors 34 mutations (including three deletions and one insertion) in the spike protein, including 15 within the RBD region (Figure 1B). These mutations are structurally focused at the top of the spike, in regions accessible to antibodies, raising the likelihood of immune evasion (Figure 1C).

### Neutralizing antibody responses to SARS-CoV-2 variants demonstrate substantial escape by Omicron

We accrued a diverse cohort of 239 COVID-19 vaccinees that were healthcare workers and/or community dwellers from Boston or Chelsea, Massachusetts (Table S1). The entire cohort had a median age of 38 years (range: 18–78 years) and was 63% female. Vaccinees had received a full series of mRNA-1273, BNT162b, or Ad26.COV2.S and were subdivided into the following subgroups: infection-naïve individuals that received their primary vaccination series within the last 3 months (“recent vax”), infection-naïve individuals that received their primary vaccination series 6–12 months before (“distant vax”), individuals that received their primary vaccination series 6–12 months before and had a history of self-reported positive PCR and/or serologic evidence of SARS-CoV-2 infection as measured by anti-nucleocapsid antibodies (“distant vax + infection”), and infection-naïve individuals that received an additional mRNA vaccine dose within the last 3 months (“booster vax”) (Figure 2A).

Sera from these cohorts were subjected to a high-throughput luminescence-based neutralization assay that we and others have previously validated (Crawford et al., 2020; Garcia-Beltran et al., 2021b; Ju et al., 2020; Moore et al., 2004; Pinto et al., 2020;

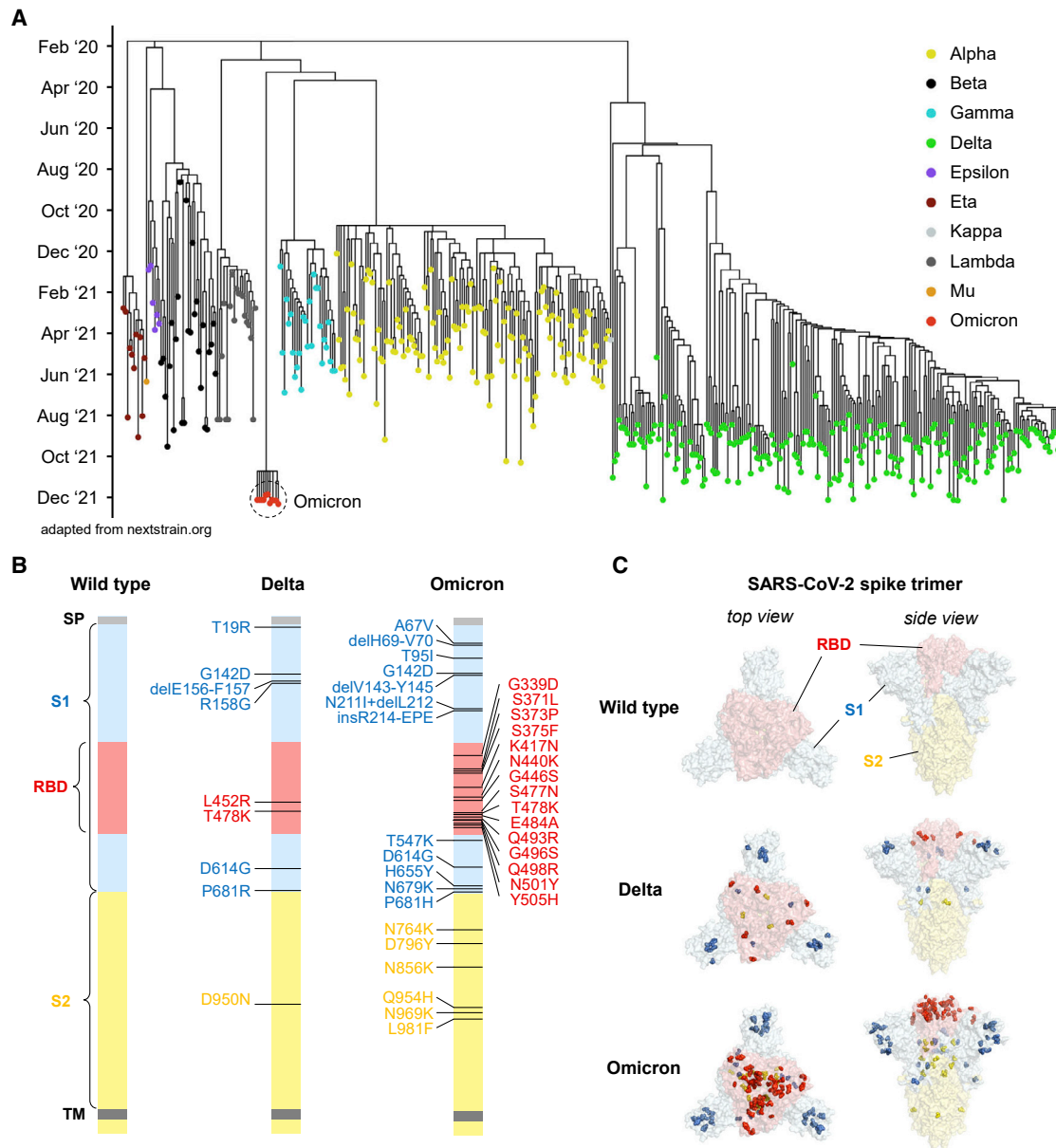

**Figure 1. Emergence of SARS-CoV-2 Omicron among global variants of concern**

(A) Phylogenetic tree of SARS-CoV-2 variants with sampling dates shows emergence of Omicron variant by December 2021 (adapted from [nextstrain.org](https://nextstrain.org); updated as of December 14, 2021).

(B) Schematic of SARS-CoV-2 spike protein structure and mutations of variants used in this study are illustrated. Omicron variant mutations used in this study were based on the most prevalent mutations (>85% frequency) found in GISAID and reflect the dominant Omicron variant. The regions within the spike protein are abbreviated as follows: SP, signal peptide; RBD, receptor binding domain; TM, transmembrane domain.


(C) Crystal structure of pre-fusion stabilized SARS-CoV-2 spike trimer (PDB ID 7JJJ) highlighting the mutational landscape of SARS-CoV-2 Delta and Omicron variants relative to SARS-CoV-2 wild type. Top views (left panels) and side views (right panels) of spike protein are shown with mutations in RBD (in red), S1 (in blue), and S2 (in yellow), highlighted with residue atoms as colored spheres.

Wang et al., 2020; Yang et al., 2020) to assess SARS-CoV-2 variant neutralization (Figure 2B). In brief, pseudovirus encoding luciferase and bearing SARS-CoV-2 variant spike proteins were exposed to dilutions of vaccinee sera prior to being added to ACE2-expressing target cells. After 48 h of co-culture, luciferase activity of the dilution series was measured to quantify infection

rates and calculate the titer that achieved 50% pseudovirus neutralization (pNT50). This pNT50 value was subsequently converted to WHO IU/mL after correcting with a WHO reference standard that was run in parallel (Knezevic et al., 2021).

In line with previous studies (Naranbhai et al., 2021a; Naranbhai et al., 2022), individuals that were recently vaccinated with



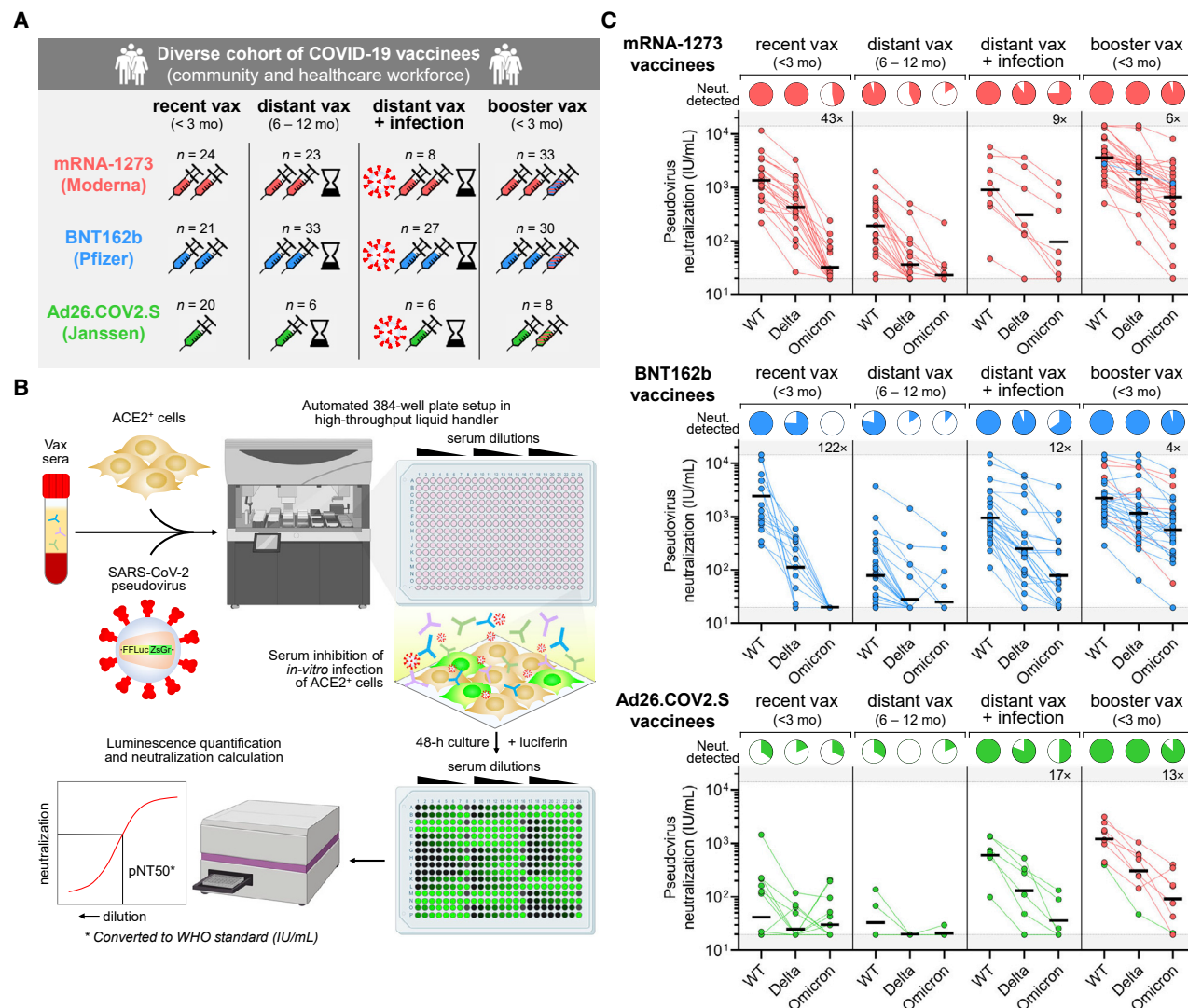

**Figure 2. Additional mRNA vaccine dose (“booster”) induces potent neutralizing responses against SARS-CoV-2 Omicron variant that are low-to-absent in primary series (“non-boostered”) vaccinees**

(A) Schematic representation of vaccinee cohorts of healthy adult community dwellers and healthcare workers. Participants who completed their primary series of vaccination with two-dose mRNA-1273 (Moderna), two-dose BNT162b (Pfizer), or 1-dose Ad26.COVS2 (Janssen) were included in this study. Vaccinees were stratified into four subgroups, as follows: infection-naïve, non-boostered individuals that received primary vaccination series within last 3 months (“recent vax”); individuals that received primary vaccination series 6–12 months before and were either without (“distant vax”) or with a history or serologic evidence of SARS-CoV-2 infection (“distant vax + infection”); and infection-naïve individuals that were boosted within the last 3 months (“booster vax”). History of SARS-CoV-2 infection was determined by either self-reported history of positive PCR test and/or positive anti-SARS-CoV-2 nucleocapsid antibody test.

(B) Schematic of experimental workflow of high-throughput SARS-CoV-2 pseudovirus neutralization assay used to determine neutralization titer of vaccinee sera against variants. Select images were created with [BioRender.com](https://www.biorender.com).

(C) Neutralization titers (in WHO IU/mL) of wild-type (WT), Delta, and Omicron pseudoviruses were determined for people who received primary vaccination series with mRNA-1273 (top panel; in red), BNT162b (middle panel; in blue), or Ad26.COVS2 (bottom panel; in green) and classified into the aforementioned subgroups (in A). Dark horizontal lines for each group denote geometric mean titer. Pie charts show the proportion of vaccinees within each group that had detectable neutralization against the indicated SARS-CoV-2 pseudovirus. Fold-decrease in geometric mean neutralization titer of Omicron relative to wild type within a subgroup is shown as a number with “x” symbol within the gray region; this was done only for vaccinee subgroups where neutralization against wild-type pseudovirus was detected in 100% of individuals. All fold-decreases shown have unadjusted  $p < 0.05$  with paired t test. Within “booster vax” subgroups (far right), boosters were homologous (same vaccine) except for 1 of 33 mRNA-1273 vaccinees that crossed-over to BNT162b (top panel; in blue), 6 of 30 BNT162b vaccinees that crossed-over to mRNA-1273 (middle panel; in red), and 7 of 8 Ad26.COVS2 vaccinees that crossed-over to mRNA-1273 (bottom panel; in red).

mRNA-1273 or BNT162b achieved substantially higher wild-type neutralization titers than Ad26.COVS2 vaccinees, with geometric mean neutralization titers (GMNTs) of 1,362 IU/mL for mRNA-

1273, 2,402 IU/mL for BNT162b, and 42 IU/mL for Ad26.COVS2 (Figure 2C). Individuals vaccinated >6 months prior exhibited substantially lower but mostly detectable wild-type neutralization

(GMNT 192 IU/mL for mRNA-1273, 78 IU/mL for BNT162b, and 33 IU/mL for Ad26.COV2.S) (Figure 2C). Prior history of infection was associated with high levels of wild-type neutralization titers even in distantly vaccinated individuals, particularly in Ad26.COV2.S vaccinees (GMNT 904 IU/mL for mRNA-1273, 947 IU/mL for BNT162b, and 603 IU/mL for Ad26.COV2.S) (Figure 2C). However, recently boosted individuals exhibited among the highest neutralization titers against wild-type SARS-CoV-2 pseudovirus (GMNT 3,862 IU/mL for mRNA-1273, 2,219 IU/mL for BNT162b, and 1,201 IU/mL for Ad26.COV2.S) (Figure 2C).

Neutralization of the Delta variant pseudovirus was decreased relative to wild type for all subgroups (Figure 2C), as has been previously reported (Mlcochova et al., 2021; Planas et al., 2021). While Delta neutralization became undetectable in most individuals vaccinated >6 months before blood draw, Delta neutralization was detectable and only modestly decreased in recently vaccinated, previously infected, and recently boosted vaccinees (Figure 2C). In contrast, Omicron neutralization was dramatically decreased among all subgroups, including recently vaccinated mRNA-1273 and BNT162b recipients, which demonstrated a complete loss of neutralization in >50% of individuals and GMNT decrease of 43-fold for mRNA-1273 and 122-fold for BNT162b (Figure 2C). Previously infected vaccinees also had a substantial decrease in Omicron neutralization titer (GMNT decrease of 9-fold for mRNA-1273, 12-fold for BNT162b, and 17-fold for Ad26.COV2.S), but most retained detectable neutralization (Figure 2C). Remarkably, however, recently boosted vaccinees exhibited potent neutralization of Omicron variant pseudovirus that was only moderately decreased relative to wild-type neutralization (GMNT decrease of 6-fold for mRNA-1273, 4-fold for BNT162b, and 13-fold for Ad26.COV2.S) (Figure 2C). Of note, among boosted vaccinees, all boosters were homologous (same vaccine) except for 1 of 33 mRNA-1273 vaccinees that crossed-over to BNT162b, 6 of 30 BNT162b vaccinees that crossed-over to mRNA-1273, and 7 of 8 Ad26.COV2.S vaccinees that crossed-over to mRNA-1273. Taken together, our results indicate that two-dose mRNA-based vaccines are effective at inducing neutralizing immunity to SARS-CoV-2 wild-type and Delta variants but suboptimal for inducing neutralizing responses to the Omicron variant.

#### Additional “booster” dose of mRNA vaccine increases breadth and cross-reactivity of neutralizing antibody response

Given the drastic increase in cross-neutralization of the SARS-CoV-2 Omicron pseudovirus in boosted versus non-boosted vaccinees, we directly compared sera from individuals that recently received their primary series to those that were boosted with an mRNA vaccine within the last 3 months. Wild-type pseudovirus neutralization was comparable between individuals who received three versus two doses of either mRNA vaccine (GMNT increase of 3-fold for mRNA-1273 and 1-fold for BNT162b). The difference in neutralization of Delta pseudoviruses between three-dose and two-dose vaccinees was similar for mRNA-1273 (GMNT increase of 3-fold) but increased in BNT162b vaccinees (GMNT increase of 9-fold). Cross-neutralization of the Omicron variant was substantially higher in individuals who received three doses of either mRNA vaccine (GMNT increase

of 19-fold for mRNA-1273 and 27-fold for BNT162b) (Figure 3A). Interestingly, Ad26.COV2.S vaccinees boosted with mRNA-1273 showed substantially higher wild-type, Delta, and Omicron pseudovirus neutralization relative to those who received Ad26.COV2.S alone (Figure 3A).


To better characterize the neutralization patterns observed between individuals who were fully vaccinated with each of the three approved vaccines and those who were boosted, we directly compared the wild-type neutralization activity of these two groups of samples against the Delta and Omicron pseudoviruses (Figure 3B). Interestingly, we found that wild-type neutralization titers from individuals who received their primary series correlated weakly to Delta variant cross-neutralization ( $R^2 = 0.35$ ; slope = 0.44;  $p < 0.0001$ ) and did not correlate with Omicron variant cross-neutralization ( $R^2 = 0.03$ ; slope = 0.05;  $p = 0.16$ ). In contrast, wild-type neutralization of boosted individuals correlated strongly with Delta ( $R^2 = 0.68$ ; slope = 1.00;  $p < 0.0001$ ) and Omicron ( $R^2 = 0.56$ ; slope = 0.94;  $p < 0.0001$ ) variant cross-neutralization. This indicates that, in addition to inducing higher neutralization titers against wild-type SARS-CoV-2, boosting increases the breadth of humoral immunity and cross-reactivity against highly mutated SARS-CoV-2 variants such as Omicron.

#### Neutralization of SARS-CoV-2 variants can be predicted by anti-spike antibody levels

Given the widespread accessibility of anti-SARS-CoV-2 spike serological assays in the clinical setting, we correlated the neutralization signal obtained by our high-throughput assay to measurements made by the Roche Elecsys anti-SARS-CoV-2 spike semi-quantitative immunoassay in the clinical laboratory. As expected, we observed higher neutralization of wild-type, Delta, and Omicron pseudoviruses with higher levels of anti-spike antibodies among all vaccinees (Figure 3C). Receiver operating characteristic (ROC) analysis was performed to assess how well anti-spike antibody levels performed at predicting positive neutralization of wild-type, Delta, and Omicron pseudoviruses, which was defined as neutralization titer >33 IU/mL (based on a previously established cohort of 1,200 pre-pandemic samples (Garcia-Beltran et al., 2021b) and converted to WHO IU/mL). Area under curve (AUC) for wild type was 0.97, Delta was 0.91, and Omicron was 0.84. The optimal cut-offs that maximized sensitivity (Se) and specificity (Sp) were determined using the “Se + Sp” method, resulting in the following: for wild type, optimal cut-off of 711 U/mL achieved 88.4% Se and 96.7% Sp; for Delta, optimal cut-off of 1,591 U/mL achieved 88.4% Se and 83.8% Sp; and for the Omicron variant, optimal cut-off of 10,300 achieved 67.2% Se and 90.6% Sp (Figures 3C and 3D). This highlights the potential use of this widely available clinical diagnostic test in predicting SARS-CoV-2 variant neutralization.

#### Omicron retains ACE2 usage and exhibits higher infection rates than other SARS-CoV-2 variants *in vitro*

Having detected a substantial degree of evasion of vaccine-induced humoral responses by Omicron, we next investigated how mutations in Omicron might result in changes in infectivity compared with wild-type SARS-CoV-2. To investigate whether the spike of circulating variants mediate host cell entry via ACE2, we examined the ability of pseudoviruses bearing wild-type, Delta,

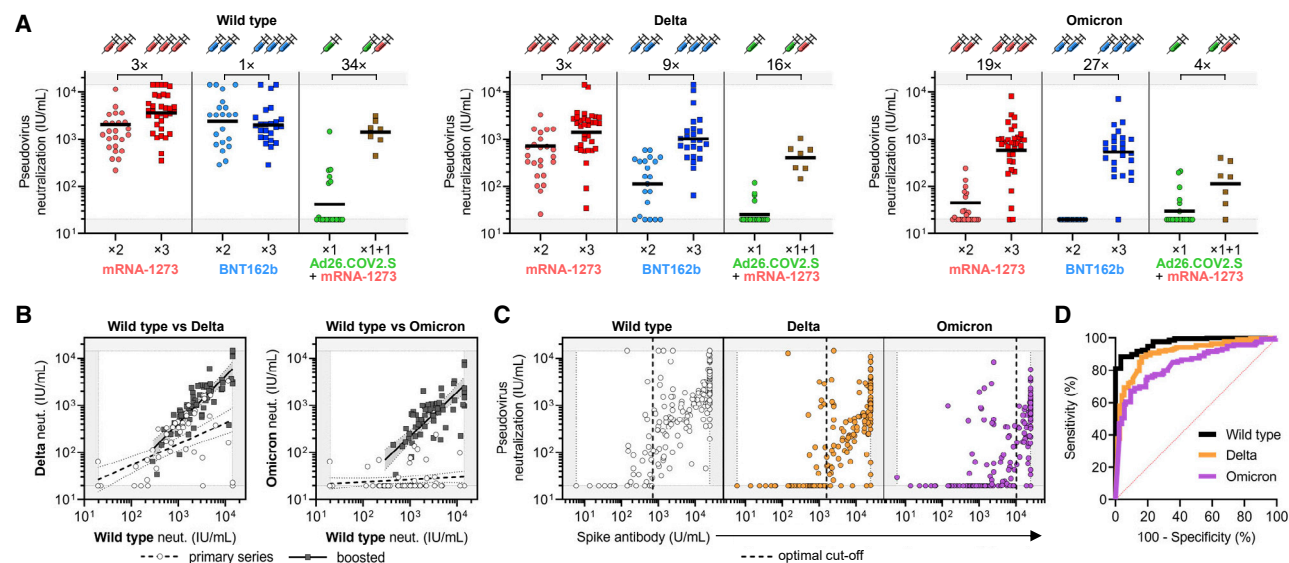

**Figure 3. Cross-reactivity of neutralizing antibody response is increased by an additional dose (“booster”) of mRNA vaccine relative to primary vaccination series and can be predicted by anti-spike antibody levels**

(A) Neutralization titers (in WHO IU/mL) of wild-type (WT; left panel), Delta (middle panel), and Omicron (right panel) SARS-CoV-2 pseudoviruses were analyzed for infection-naïve participants that were recently vaccinated with primary series or booster (<3 months). Recently vaccinated individuals received mRNA-1273 (x2), BNT162b (x2), or Ad26.COV2.S (x1), and boosted individuals received a homologous booster of mRNA-1273 (x3), or BNT162b (x3), or a cross-over booster of mRNA-1273 for Ad26.COV2.S vaccinees (x1 + 1). Fold-increase in geometric mean neutralization titer of boosted versus non-boosted individuals is shown as a number with “x” symbol. This analysis is based on experimental data depicted in Figure 2C but excludes participants with prior SARS-CoV-2 infection, distant vaccination (>6 months), and/or cross-over between mRNA-1273 and BNT162b to understand differences in neutralizing responses soon after primary vaccination series versus boosting. The single Ad26.COV2.S vaccinee that received a homologous booster with Ad26.COV2.S was also excluded from this analysis. Dark horizontal lines for each group denote geometric mean titer.

(B) Aggregate data from study participants (A) that recently received primary vaccination series (“primary series”; white circles) or were recently boosted (“boosted”; dark gray squares) were used for linear regression analysis of wild type versus Delta (left panel) or wild type versus Omicron (right panel) pseudovirus neutralization. Wild-type neutralization titers correlated with Delta neutralization in “primary series” individuals ( $R^2 = 0.35$ ; slope = 0.44;  $p < 0.0001$ ) and even more strongly in “boosted” individuals ( $R^2 = 0.68$ ; slope = 1.00;  $p < 0.0001$ ). Wild-type neutralization titers showed no significant relationship with Omicron neutralization in “primary series” individuals ( $R^2 = 0.03$ ; slope = 0.05;  $p = 0.16$ ); however, “boosted” individuals showed a significant correlation with Omicron neutralization titers ( $R^2 = 0.56$ ; slope = 0.94;  $p < 0.0001$ ).

(C) Anti-SARS-CoV-2 spike antibodies levels (measured by an EUA-approved clinical diagnostic test) of all vaccinees were plotted against neutralization of wild-type (left panel; white circles), Delta (middle panel; orange circles), and Omicron (right panel; purple circles) SARS-CoV-2 pseudoviruses. Optimal spike antibody cut-offs for predicting positive neutralization were determined by ROC analyses in (D) and are indicated with a vertical dashed line.

(D) Receiver operating characteristic (ROC) analyses assessing the ability of spike antibody levels to predict neutralization of wild-type (black line), Delta (orange line), and Omicron (purple line) pseudoviruses. Positive neutralization was defined as >33 IU/mL (previously defined with a cohort of 1,200 pre-pandemic samples (Garcia-Beltran et al., 2021b) and converted to WHO IU/mL). Area under curve (AUC) for wild type was 0.97, Delta was 0.91, and Omicron was 0.84, with  $p < 0.0001$  for all three. Optimal cut-offs that maximized sensitivity (Se) and specificity (Sp) were determined using the “Se + Sp” method, and were as follows: for wild type, optimal cut-off of 711 U/mL achieved 88.4% Se and 96.7% Sp; for Delta, optimal cut-off of 1,591 U/mL achieved 88.4% Se and 83.8% Sp; and for Omicron variant, optimal cut-off of 10,300 achieved 67.2% Se and 90.6% Sp. These optimal cut-off values are plotted as a vertical dashed line in (C).


and Omicron spike to infect 293T-ACE2 cells or parental 293T cells devoid of ACE2 receptor (Figure 4A). None of the pseudoviruses tested infected 293T cells in the absence of ACE2 protein, confirming that, like other SARS-CoV-2 strains, Omicron uses ACE2 for host cell entry (Figure 4B). To comprehensively determine quantitative differences in the infectivity of SARS-CoV-2 variants of concern, we compared the efficiency of a panel of pseudovirus variants to infect 293T-ACE2 cells over a range of viral concentrations. Remarkably, Omicron pseudovirus exhibited greater infection of target cells regardless of concentration when compared with all other tested variants (Figure 4C). Comparison of the linear regressions of each pseudovirus to wild type over the entire range revealed that, whereas the Gamma variant exhibited similar infection rates to wild type, Beta showed less infec-

tion, and Delta was nearly 2-fold more efficient at infecting target cells. Strikingly, Omicron exhibited infection rates that were 4-fold higher than wild type and 2-fold higher than Delta. Taken together, these data strongly suggest distinct differences in infectivity according to spike sequence, with Omicron exhibiting more efficient ACE2-mediated infection than wild type or other variant strains (Figure 4D).

## DISCUSSION

Prior characterization of vaccine-induced humoral immunity against SARS-CoV-2 variants of concern revealed significant loss of activity against the Gamma and Beta variants, owing largely to three mutations in the RBD region of the spike



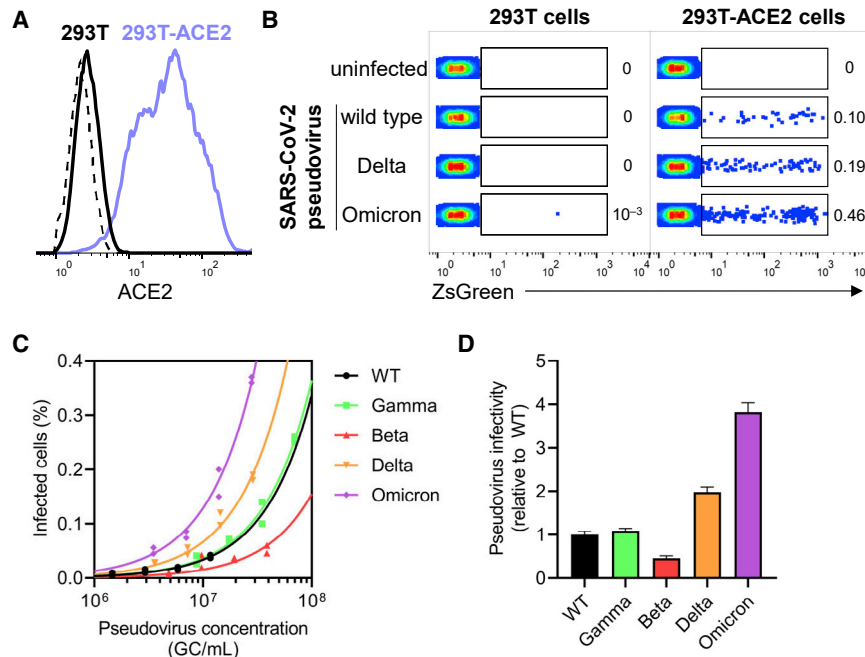

**Figure 4. SARS-CoV-2 Omicron pseudovirus demonstrates a substantial increase in infectivity of ACE2<sup>+</sup> cells relative to other SARS-CoV-2 variants *in vitro***

(A) Flow cytometry histogram depicting ACE2 surface staining on parental 293T cell line (in black), and 293T-ACE2 cell line stably expressing human ACE2 (in lilac). Dashed line indicates unstained control.

(B) Representative dot plots show percentage of 293T (left panel) and 293T-ACE2 (right panel) cells that were infected with wild-type, Delta, or Omicron SARS-CoV-2 pseudoviruses (within ZsGreen<sup>+</sup> gate) after 48 h of co-culture. Pseudoviruses were produced in parallel and under identical conditions.

(C) Titration of SARS-CoV-2 pseudoviruses of wild-type, Gamma, Beta, Delta, and Omicron variants was performed on 293T-ACE2 cells and correlated to pseudovirus concentration in genome copies (GC)/mL as determined by qPCR. Two technical replicates ( $n = 2$ ) and a linear regression were performed to fit data, given that a linear relationship of virus concentration versus infected cells can be assumed at infection rates <10%.

(D) Pseudovirus infectivity relative to wild type was measured for each SARS-CoV-2 variant in (C) by calculating fold change in slope from (C) for each pseudovirus relative to wild type. Bars and error bars depict mean and standard error of the mean.

(Garcia-Beltran et al., 2021a). Given the 15 mutations present in the RBD of Omicron, which overlap the three sites mutated in Gamma and Beta (K417, E484, and N501), as well as prior work done to understand the potential for antibody escape (Schmidt et al., 2021), it was anticipated that this variant would be significantly more resistant to neutralization by vaccinee sera. To experimentally validate this assumption, we generated a pseudovirus with the Omicron spike incorporating 34 distinct mutations (including three deletions and one insertion) relative to the Wuhan spike. Importantly, our spike did not contain R346K, which has been observed at relatively low frequency (<10%) within the Omicron lineage, but which is known to mediate escape from class 3 neutralizing monoclonal antibodies such as AZD1061/cilgavimab (Barnes et al., 2020; Dong et al., 2021; Greaney et al., 2021b).

As stated in our previous study, given the assumption of a polyclonal response to vaccination, we would have initially anticipated that small numbers of variations in spike protein would have only modest effects on recognition by the immune system (Garcia-Beltran et al., 2021a). However, our earlier findings suggested that the response to primary vaccination was insufficient to effectively neutralize spikes containing a relatively small number of mutations. Consistent with this observation, we found that the Omicron variant, harboring substantially more mutations than prior variants, efficiently escapes humoral immunity induced by primary vaccination. Strikingly, however, additional mRNA vaccine doses appear to enable cross-neutralizing responses against Omicron, either by further affinity maturation of existing antibodies or targeting of new epitopes shared among variants.

Our study included distantly vaccinated individuals (>6 months) with evidence of prior infection. However, our approach could not distinguish between infection prior to vaccination or asymptomatic breakthrough infection. Irrespective of this distinction, these indi-

viduals exhibited higher levels of neutralizing activity against the wild-type, Delta, and Omicron pseudoviruses as compared with distantly vaccinated, infection-naïve individuals. These results suggest that additional exposure to viral antigens through infection yields enhanced overall neutralization activity. It is important to note, however, that an additional dose of mRNA vaccine in infection-naïve vaccinees yielded substantially higher cross-neutralizing activity against Omicron as compared with prior infection.

Interestingly, vaccinated individuals who received an additional “booster” mRNA vaccine dose generated the most potent and cross-reactive antibody responses of any cohort tested, including individuals who recently received their primary series. The effect was similar to that observed following mRNA vaccine boosting in rhesus macaques (Corbett et al., 2021). This substantial difference in Omicron cross-neutralization was surprising, given that two-dose and three-dose vaccine regimens elicited comparable levels of wild-type neutralization and because the immunogen was identical for all doses, which would be expected to elicit memory responses from the existing, strain-specific antibody repertoire. This difference in cross-reactivity between vaccinees that recently received three versus two doses of vaccine could be the consequence of a third dose amplifying existing, low-level neutralizing antibodies that target conserved epitopes on spike protein or increasing the affinity of existing neutralizing antibodies so as to render them less sensitive to mutations in their target epitopes, or both. Although the emergence of variants has catalyzed the development of variant-specific “booster” shots to increase variant neutralization, our results suggest that current wild-type-based mRNA vaccines are sufficient to stimulate cross-reactive humoral responses greater than might have been anticipated. Whether this is the result of third-dose timing or because the primary series was



simply insufficient to fully stimulate broadly neutralizing B cell responses to spike antigen is unclear.

We also found marked differences between approved vaccine regimens, with mRNA-1273 and BNT162b2 vaccines offering significantly greater humoral immune responses against all variants than Ad26.COV2.S, as has been previously described (Naranbhai et al., 2021a; Naranbhai et al., 2022). Recipients of mRNA vaccines following Ad26.COV2.S exhibited markedly improved neutralizing titers but did not achieve the breadth of response seen with individuals that received three doses of mRNA vaccine. Our results would suggest that these recipients of Ad26.COV2.S vaccines may benefit from additional mRNA vaccine doses with the potential to further raise titers and broaden their neutralizing activity.

Several reports have measured the decay of neutralizing antibody responses raised by mRNA COVID-19 vaccines with a half-life of 69–173 days (Doria-Rose et al., 2021; Levin et al., 2021). Our results from distantly vaccinated samples are in line with these reports, given that neutralization activity against wild-type pseudovirus was more than 10-fold lower in distantly vaccinated individuals. Given the activity generated against Omicron, it will be important to determine the longevity of cross-neutralizing antibody responses in individuals that received a third dose of mRNA vaccine, as this could have a meaningful impact on vaccine efficacy as future variants of concern continue to emerge.

In support of the rapid spread of Omicron globally (BBC News, 2021), our *in vitro* studies strongly suggest that Omicron spike mediates more efficient target cell entry, resulting in increased infectivity. Interestingly, while the Beta variant previously exhibited escape from vaccine-induced humoral immunity, we found that infectivity of Beta was lower than wild type, perhaps explaining the relatively lower epidemic spread of Beta. In contrast, the now globally dominant Delta strain exhibited 2-fold increased efficiency of pseudovirus entry relative to wild type, consistent with prior studies reporting more efficient spike fusion that may have contributed to its rapid spread (Zhang et al., 2021). Additional studies examining the interaction between Omicron spike and ACE2, as well as other potential interaction partners, will be necessary to discern whether more efficient infection of the SARS-CoV-2 pseudovirus is also observed in replication-competent viruses.

Taken together, we demonstrate that Omicron drastically escapes vaccine-induced immunity after primary vaccination series with mRNA-1273 (Moderna), BNT162b2 (Pfizer-BioNTech), or Ad26.COV2.S (Johnson & Johnson/Janssen) and exhibits increased pseudovirus infection rates *in vitro*, raising the potential for increased transmissibility. It is worth noting that, despite escape from humoral immunity, Omicron breakthrough infections may result in attenuated disease severity in vaccinees due to pre-existing cellular and innate immunity. However, neutralization remains the leading correlate of protection from infection, and this study demonstrates that receiving a third dose of an mRNA-based vaccine effectively yields a potent cross-neutralizing response against SARS-CoV-2 Omicron, likely through increasing breadth and cross-reactivity of neutralizing antibodies. These findings support the need for rapid and synchronized widespread deployment of additional mRNA vac-

cine doses as a public health measure to curtail the emergence and spread of highly mutated SARS-CoV-2 variants.

### Limitations of the study

Although previous studies that used pseudovirus neutralization to model the sensitivity of replicating SARS-CoV-2 to neutralizing antibodies have shown excellent correlations (Crawford et al., 2020; Ju et al., 2020; Moore et al., 2004; Pinto et al., 2020; Riepler et al., 2020; Wang et al., 2020; Yang et al., 2020), it is possible that the mutations in the Omicron spike protein may cause the Omicron pseudovirus to behave differently than previously tested variants. However, recent reports have demonstrated similar loss of neutralizing activity by vaccinee sera against intact Omicron coronavirus (Cele et al., 2021b). In addition, while we confirmed that ACE2 expression is required for infection of 293T cells, natural target cells in the respiratory tract may express alternative receptors or attachment factors that facilitate infection and are not adequately modeled in our system. In addition, our cohort was cross-sectional and not longitudinal, which limits our ability to estimate changes in neutralization titers over time across single individuals. Furthermore, we did not assess other antibody-mediated functions, such as complement deposition, antibody-dependent cellular cytotoxicity, or antibody-dependent cellular phagocytosis, which may contribute to protection even in the absence of neutralizing antibodies. We did not assess the role of vaccine-elicited cellular immune responses mediated by T cells and NK cells, which are likely to play a key role in disease prevention for vaccine recipients.

### STAR★METHODS

Detailed methods are provided in the online version of this paper and include the following:

- KEY RESOURCES TABLE
- RESOURCE AVAILABILITY
  - Lead contact
  - Materials availability
  - Data and code availability
- EXPERIMENTAL MODEL AND SUBJECT DETAILS
  - Human subjects
  - Cell lines
- METHOD DETAILS
  - Construction of variant spike expression plasmids
  - SARS-CoV-2 pseudovirus neutralization assay
  - Titering and flow cytometry
  - Quantitation of pseudovirus by RT-qPCR
  - Anti-SARS-CoV-2 spike and nucleocapsid antibody assays
- QUANTIFICATION AND STATISTICAL ANALYSIS

### SUPPLEMENTAL INFORMATION

Supplemental information can be found online at <https://doi.org/10.1016/j.cell.2021.12.033>.

## ACKNOWLEDGMENTS

We thank Anand Dighe, MD; Andrea Nixon, BS; and the MGH Core Laboratory for excellent assistance with clinical SARS-CoV-2 serology testing. We wish to thank Michael Farzan, PhD, for providing ACE2-expressing 293T cells. A.H. is supported by the DZIF (German Center for Infection Research, TTU 01.709; 8009701709). B.M.H. is supported by award number T32GM007753 from the National Institute of General Medical Sciences and F30 AI160908 from the National Institute of Allergy and Infectious Diseases. J.F. is supported by T32AI007245. D.J.G. and M.C.P. were supported by the VIC Innovation Fund. A.G.S. was supported by NIH R01 AI146779 and a Massachusetts Consortium on Pathogenesis Readiness (MassCPR) grant. A.H. is supported by the DZIF (German Center for Infection Research, TTU 01.709). This work was supported by the Peter and Ann Lambertus Family Foundation. V.N. received support from a Medscape Young Investigators Lung Cancer Award. A.B.B. was supported by the National Institute on Drug Abuse (NIDA) Avenir New Innovator Award DP2DA040254, the MGH Transformative Scholars Program, and a Massachusetts Consortium on Pathogenesis Readiness (MassCPR) grant.

## AUTHOR CONTRIBUTIONS

W.F.G.B., K.J.S.D., E.C.L., V.N., and A.B.B. designed the experiments. W.F.G.B., K.J.S.D., E.C.L., M.L.S., A.D.N., and A.B.B. carried out experiments and analyzed data. W.F.G.B., V.N., C.B., O.O., C.C.C., D.J.G., M.C.P., and A.J.I. led and performed cohort recruitment and sampling. J.F., B.M.H., and A.G.S. provided key reagents and useful discussions and insights. A.L.R. performed sequence validation of plasmids. A.B.B., V.N., W.F.G.B., A.D.N., and A.H. contributed to statistical and sequence analyses. A.J.I. and V.N. provided key discussions and input into experimental design. W.F.G.B., V.N., A.H., and A.B.B. wrote the paper with contributions from all authors.

## DECLARATION OF INTERESTS

The authors declare no competing interests.

## INCLUSION AND DIVERSITY

We worked to ensure gender balance in the recruitment of human subjects. We worked to ensure ethnic or other types of diversity in the recruitment of human subjects. We worked to ensure that the study questionnaires were prepared in an inclusive way. One or more of the authors of this paper self-identifies as an underrepresented ethnic minority in science. One or more of the authors of this paper self-identifies as a member of the LGBTQ+ community.

Received: December 14, 2021

Revised: December 18, 2021

Accepted: December 20, 2021

Published: December 23, 2021

## REFERENCES

Baden, L.R., El Sahly, H.M., Essink, B., Kottloff, K., Frey, S., Novak, R., Diemert, D., Spector, S.A., Roupheal, N., Creech, C.B., et al. (2021). Efficacy and safety of the mRNA-1273 SARS-CoV-2 vaccine. *N. Engl. J. Med.* **384**, 403–416.

Bajema, K.L., Dahl, R.M., Evener, S.L., Prill, M.M., Rodriguez-Barradas, M.C., Marconi, V.C., Beenhouwer, D.O., Holodniy, M., Lucero-Obusan, C., Brown, S.T., et al. (2021). Comparative effectiveness and antibody responses to Moderna and Pfizer-BioNTech COVID-19 vaccines among hospitalized veterans—five Veterans Affairs Medical Centers, United States, February 1–September 30, 2021. *MMWR Morb. Mortal. Wkly. Rep.* **70**, 1700–1705.

Barnes, C.O., Jette, C.A., Abernathy, M.E., Dam, K.-M.A., Esswein, S.R., Grinstead, H.B., Malyutin, A.G., Sharaf, N.G., Huey-Tubman, K.E., Lee, Y.E., et al. (2020). SARS-CoV-2 neutralizing antibody structures inform therapeutic strategies. *Nature* **588**, 682–687.

Bar-On, Y.M., Goldberg, Y., Mandel, M., Bodenheimer, O., Freedman, L., Kalkstein, N., Mizrahi, B., Alroy-Preis, S., Ash, N., Milo, R., and Huppert, A. (2021). Protection of BNT162b2 vaccine booster against Covid-19 in Israel. *N. Engl. J. Med.* **385**, 1393–1400.

Cele, S., Gazy, I., Jackson, L., Hwa, S.-H., Tegally, H., Lustig, G., Giandhari, J., Pillay, S., Wilkinson, E., Naidoo, Y., et al. (2021a). Escape of SARS-CoV-2 501Y.V2 from neutralization by convalescent plasma. *Nature* **593**, 142–146.

Cele, S., Jackson, L., Khan, K., Khoury, D.S., Moyo-Gwete, T., Tegally, H., Scheepers, C., Amoako, D., Karim, F., Bernstein, M., et al. (2021b). SARS-CoV-2 Omicron has extensive but incomplete escape of Pfizer BNT162b2 elicited neutralization and requires ACE2 for infection. *medRxiv*. <https://doi.org/10.1101/2021.12.08.21267417>.

BBC News (2021). Covid: Omicron spreading at unprecedented rate, WHO says. <https://www.bbc.com/news/world-59656385>.

Corbett, K.S., Gagne, M., Wagner, D.A., O'Connell, S., Narpala, S.R., Flebbe, D.R., Andrew, S.F., Davis, R.L., Flynn, B., Johnston, T.S., et al. (2021b). Protection against SARS-CoV-2 Beta variant in mRNA-1273 vaccine–boosted nonhuman primates. *Science* **374**, 6573. <https://doi.org/10.1126/science.abl8912>.

Crawford, K.H.D., Eguia, R., Dingens, A.S., Loes, A.N., Malone, K.D., Wolf, C.R., Chu, H.Y., Tortorici, M.A., Veesler, D., Murphy, M., et al. (2020). Protocol and reagents for pseudotyping lentiviral particles with SARS-CoV-2 spike protein for neutralization assays. *Viruses* **12**.

Cromer, D., Steain, M., Reynaldi, A., Schlub, T.E., Wheatley, A.K., Juno, J.A., Kent, S.J., Triccas, J.A., Khoury, D.S., and Davenport, M.P. (2021). Neutralizing antibody titres as predictors of protection against SARS-CoV-2 variants and the impact of boosting: a meta-analysis. *The Lancet Microbe*. [https://doi.org/10.1016/S2666-5247\(21\)00267-6](https://doi.org/10.1016/S2666-5247(21)00267-6).

Dong, J., Zost, S.J., Greaney, A.J., Starr, T.N., Dingens, A.S., Chen, E.C., Chen, R.E., Case, J.B., Sutton, R.E., Gilchuk, P., et al. (2021). Genetic and structural basis for SARS-CoV-2 variant neutralization by a two-antibody cocktail. *Nat. Microbiol.* **6**, 1233–1244.

Doria-Rose, N., Suthar, M.S., Makowski, M., O'Connell, S., McDermott, A.B., Flach, B., Ledgerwood, J.E., Mascola, J.R., Graham, B.S., Lin, B.C., et al. (2021). Antibody persistence through 6 months after the second dose of mRNA-1273 vaccine for Covid-19. *N. Engl. J. Med.* **384**, 2259–2261.

Earle, K.A., Ambrosino, D.M., Fiore-Gartland, A., Goldblatt, D., Gilbert, P.B., Siber, G.R., Dull, P., and Plotkin, S.A. (2021). Evidence for antibody as a protective correlate for COVID-19 vaccines. *Vaccine* **39**, 4423–4428.

Garcia-Beltran, W.F., Lam, E.C., Astudillo, M.G., Yang, D., Miller, T.E., Feldman, J., Hauser, B.M., Caradonna, T.M., Clayton, K.L., Nitido, A.D., et al. (2021b). COVID-19-neutralizing antibodies predict disease severity and survival. *Cell* **184**, 476–488.e11.

Garcia-Beltran, W.F., Lam, E.C., St. Denis, K., Nitido, A.D., Garcia, Z.H., Hauser, B.M., Feldman, J., Pavlovic, M.N., Gregory, D.J., Poznansky, M.C., et al. (2021a). Multiple SARS-CoV-2 variants escape neutralization by vaccine-induced humoral immunity. *Cell* **184**, 2372–2383.e9.

Greaney, A.J., Loes, A.N., Crawford, K.H.D., Starr, T.N., Malone, K.D., Chu, H.Y., and Bloom, J.D. (2021a). Comprehensive mapping of mutations in the SARS-CoV-2 receptor-binding domain that affect recognition by polyclonal human plasma antibodies. *Cell Host Microbe* **29**, 463–476.e6.

Greaney, A.J., Starr, T.N., Barnes, C.O., Weisblum, Y., Schmidt, F., Caskey, M., Gaebler, C., Cho, A., Agudelo, M., Finkin, S., et al. (2021b). Mapping mutations to the SARS-CoV-2 RBD that escape binding by different classes of antibodies. *Nat. Commun.* **12**, 4196.

Ju, B., Zhang, Q., Ge, J., Wang, R., Sun, J., Ge, X., Yu, J., Shan, S., Zhou, B., Song, S., et al. (2020). Human neutralizing antibodies elicited by SARS-CoV-2 infection. *Nature* **584**, 115–119.

Khoury, D.S., Cromer, D., Reynaldi, A., Schlub, T.E., Wheatley, A.K., Juno, J.A., Subbarao, K., Kent, S.J., Triccas, J.A., and Davenport, M.P. (2021). Neutralizing antibody levels are highly predictive of immune protection from symptomatic SARS-CoV-2 infection. *Nat. Med.* **27**, 1205–1211.

- Knezevic, I., Mattiuzzo, G., Page, M., Minor, P., Griffiths, E., Nuebling, M., and Moorthy, V. (2021). WHO International Standard for evaluation of the antibody response to COVID-19 vaccines: call for urgent action by the scientific community. *Lancet Microbe*. [https://doi.org/10.1016/S2666-5247\(21\)00266-4](https://doi.org/10.1016/S2666-5247(21)00266-4).
- Levin, E.G., Lustig, Y., Cohen, C., Fluss, R., Indenbaum, V., Amit, S., Doolman, R., Asraf, K., Mendelson, E., Ziv, A., et al. (2021). Waning immune humoral response to BNT162b2 Covid-19 vaccine over 6 months. *N. Engl. J. Med.* **385**, e84.
- Lopez Bernal, J., Andrews, N., Gower, C., Gallagher, E., Simmons, R., Thelwall, S., Stowe, J., Tessier, E., Groves, N., Dabrera, G., et al. (2021). Effectiveness of Covid-19 vaccines against the B.1.617.2 (Delta) variant. *N. Engl. J. Med.* **385**, 585–594.
- Micochova, P., Kemp, S.A., Dhar, M.S., Papa, G., Meng, B., Ferreira, I.A.T.M., Datir, R., Collier, D.A., Albecka, A., Singh, S., et al. (2021). SARS-CoV-2 B.1.617.2 Delta variant replication and immune evasion. *Nature* **599**, 114–119.
- Moore, M.J., Dorfman, T., Li, W., Wong, S.K., Li, Y., Kuhn, J.H., Coderre, J., Vasilieva, N., Han, Z., Greenough, T.C., et al. (2004). Retroviruses pseudotyped with the severe acute respiratory syndrome coronavirus spike protein efficiently infect cells expressing angiotensin-converting enzyme 2. *J. Virol* **78**, 10628–10635.
- Naranbhai, V., Garcia-Beltran, W.F., Chang, C.C., Mairena, C.B., Thierauf, J.C., Kirkpatrick, G., Onozato, M.L., Cheng, J., St Denis, K.J., Lam, E.C., et al. (2021a). Comparative immunogenicity and effectiveness of mRNA-1273, BNT162b2 and Ad26.COV2.S COVID-19 vaccines. *J. Infect. Dis.*
- Naranbhai, V., St. Denis, K.J., Lam, E.C., Ofoman, O., Garcia-Beltran, W.F., Berrios, C., Bhan, A.K., Gainor, J.F., Balazs, A.B., and Iafrate, A.J. (2021b). Neutralization breadth of SARS-CoV-2 viral variants following primary series and booster SARS-CoV-2 vaccines in patients with cancer. *Cancer Cell*.
- Naranbhai, V., Pernat, C.A., Gavralidis, A., St Denis, K.J., Lam, E.C., Spring, L.M., Isakoff, S.J., Farmer, J.R., Zubiri, L., Hobbs, G.S., et al. (2022). Immunogenicity and reactivity of SARS-CoV-2 vaccines in patients with cancer: the CANVAX cohort study. *J. Clin. Oncol.* **40**, 12–23.
- Pinto, D., Park, Y.-J., Beltramello, M., Walls, A.C., Tortorici, M.A., Bianchi, S., Jaconi, S., Culp, K., Zatta, F., De Marco, A., et al. (2020). Structural and functional analysis of a potent sarbecovirus neutralizing antibody. *bioRxiv*. <https://doi.org/10.1101/2020.04.07.023903>.
- Planas, D., Veyer, D., Baidaliuk, A., Staropoli, I., Guivel-Benhassine, F., Rajah, M.M., Planchais, C., Porrot, F., Robillard, N., Puech, J., et al. (2021). Reduced sensitivity of SARS-CoV-2 variant Delta to antibody neutralization. *Nature* **596**, 276–280.
- Polack, F.P., Thomas, S.J., Kitchin, N., Absalon, J., Gurtman, A., Lockhart, S., Perez, J.L., Pérez Marc, G., Moreira, E.D., Zerbini, C., et al. (2020). Safety and efficacy of the BNT162b2 mRNA Covid-19 vaccine. *N. Engl. J. Med.* **383**, 2603–2615.
- Rambaut, A., Holmes, E.C., O’Toole, Á., Hill, V., McCrone, J.T., Ruis, C., du Plessis, L., and Pybus, O.G. (2020). A dynamic nomenclature proposal for SARS-CoV-2 lineages to assist genomic epidemiology. *Nat. Microbiol.* **5**, 1403–1407.
- Riepler, L., Rössler, A., Falch, A., Volland, A., Borena, W., von Laer, D., and Kimpel, J. (2020). Comparison of four SARS-CoV-2 neutralization assays. *Vaccines (Basel)* **9**, 13.
- Sadoff, J., Gray, G., Vandebosch, A., Cárdenas, V., Shukarev, G., Grinsztejn, B., Goepfert, P.A., Truyers, C., Fennema, H., Spiessens, B., et al. (2021). Safety and efficacy of single-dose Ad26.COV2.S vaccine against Covid-19. *N. Engl. J. Med.* **384**, 2187–2201.
- Schmidt, F., Weisblum, Y., Rutkowska, M., Poston, D., DaSilva, J., Zhang, F., Bednarski, E., Cho, A., Schaefer-Babajew, D.J., Gaebler, C., et al. (2021). High genetic barrier to SARS-CoV-2 polyclonal neutralizing antibody escape. *Nature* **600**, 512–516.
- Siebring-van Olst, E., Vermeulen, C., de Menezes, R.X., Howell, M., Smit, E.F., and van Beusechem, V.W. (2013). Affordable luciferase reporter assay for cell-based high-throughput screening. *J. Biomol. Screen* **18**, 453–461.
- Tian, F., Tong, B., Sun, L., Shi, S., Zheng, B., Wang, Z., Dong, X., and Zheng, P. (2021). N501Y mutation of spike protein in SARS-CoV-2 strengthens its binding to receptor ACE2. *eLife* **10**, e69091.
- Tregoning, J.S., Flight, K.E., Higham, S.L., Wang, Z., and Pierce, B.F. (2021). Progress of the COVID-19 vaccine effort: viruses, vaccines and variants versus efficacy, effectiveness and escape. *Nat. Rev. Immunol.* **21**, 626–636.
- Wang, C., Li, W., Drabek, D., Okba, N.M.A., van Haperen, R., Osterhaus, A.D.M.E., van Kuppeveld, F.J.M., Haagmans, B.L., Grosveld, F., and Bosch, B.-J. (2020). A human monoclonal antibody blocking SARS-CoV-2 infection. *Nat. Commun.* **11**, 2251.
- Yang, R., Huang, B., Ruhan, A., Li, W., Wang, W., Deng, Y., and Tan, W. (2020). Development and effectiveness of Pseudotyped SARS-CoV-2 system as determined by neutralizing efficiency and entry inhibition test in vitro. *Biosaf. Health* **2**, 226–231.
- Zhang, J., Xiao, T., Cai, Y., Lavine, C.L., Peng, H., Zhu, H., Anand, K., Tong, P., Gautam, A., Mayer, M.L., et al. (2021). Membrane fusion and immune evasion by the spike protein of SARS-CoV-2 Delta variant. *Science* **374**, 1353–1360.
- Zhou, D., Dejnirattisai, W., Supasa, P., Liu, C., Mentzer, A.J., Ginn, H.M., Zhao, Y., Duyvesteyn, H.M.E., Tuekprakhon, A., Nutalai, R., et al. (2021). Evidence of escape of SARS-CoV-2 variant B.1.351 from natural and vaccine-induced sera. *Cell* **184**, 2348–2361.e6.

STAR★METHODS

KEY RESOURCES TABLE

REAGENT or RESOURCE	SOURCE	IDENTIFIER
<b>Antibodies</b>		
Human ACE-2 Alexa Fluor 647-conjugated Antibody	R&D	FAB9332R-100G; Lot #: 1580709
<b>Bacterial and virus strains</b>		
SARS-CoV-2 pseudotyped pHAGE-CMV-Luc2-IRES-ZsGreen-W lentivirus	This study	N/A
SARS-CoV-2 ΔC18 (wild type) pseudotyped pHAGE-CMV-Luc2-IRES-ZsGreen-W lentivirus	This study	N/A
SARS-CoV-2 ΔC18 B.1.617.2 pseudotyped pHAGE-CMV-Luc2-IRES-ZsGreen-W lentivirus	This study	N/A
SARS-CoV-2 ΔC18 B.1.351 pseudotyped pHAGE-CMV-Luc2-IRES-ZsGreen-W lentivirus	This study	N/A
SARS-CoV-2 ΔC18 P.1 pseudotyped pHAGE-CMV-Luc2-IRES-ZsGreen-W lentivirus	This study	N/A
DH5α Zymo-Competent <i>E. Coli</i>	Zymo	Cat# T3009
<b>Chemicals, peptides, and recombinant proteins</b>		
Phosphate buffered saline (PBS)	Corning	Cat# 21-031-CV
Dulbecco's modified eagle medium (DMEM)	Corning	Cat# 10-013-CV
Fetal bovine serum (FBS)	VWR	Cat# 89510-186
Penicillin/streptomycin	Corning	Cat# 30-002-CI
Polyethylenimine 25K MW, linear	Polysciences Inc.	Cat# 23966
Puromycin	Sigma	Cat# P8833-10MG
ATP	Sigma	Cat# A2383-5G
Magnesium chloride	BDH	Cat# BDH9244-500G
Magnesium sulfate	BDH	Cat# BDH9246-500G
Dithiothreitol (DTT)	VWR	Cat# 97061-338
D-luciferin	Gold Bio	Cat# LUCK-2G
EDTA	Sigma	Cat# 03690-100ML
Triton-X100	Fisher	Cat# BP151-500
Polybrene	Sigma	Cat# H9268-5G
qScript XLT 1-Step RT-qPCR ToughMix, Low-ROX	Quantabio	Cat# 95134-500
Turbo Dnase	Fisher	Cat# AM2239
<b>Experimental models: Cell lines</b>		
HEK293T/17 Cells	ATCC	Cat# CRL-11268
293T/ACE2.MF	Obtained from the lab of Dr. Michael Farzan	N/A

(Continued on next page)



**Continued**

REAGENT or RESOURCE	SOURCE	IDENTIFIER
<b>Oligonucleotides</b>		
Lentivirus LTR qPCR Probe /56-FAM/AGC+C/i5NitInd/GG+GA/ZEN/GCTCTCTGGC/3IABkFQ/	Integrated DNA Technology	N/A
Lentivirus LTR qPCR 5' Primer GGTCTCTCTIGITAGACCAG	Integrated DNA Technology	N/A
Lentivirus LTR qPCR 3' Primer TTTATTGAGGCTTAAGCAGTGGG	Integrated DNA Technology	N/A
<b>Recombinant DNA</b>		
pHAGE-CMV-Luc2-IRES-ZsGreen-W (backbone)	This study	Addgene Cat# 164432
pTwist-SARS-CoV-2 Δ18 (spike)	This study	Addgene Cat# 164436
pTwist-SARS-CoV-2 Δ18 B.1.617.2v1 (spike)	This study	Addgene Cat# 179905
pTwist-SARS-CoV-2 Δ18 B.1.351v2 (spike)	This study	Addgene Cat# 169463
pTwist-SARS-CoV-2 Δ18 P.1 (spike)	This study	Addgene Cat# 173476
pTwist-SARS-CoV-2 Δ18 B.1.1.529 (spike)	This study	Addgene Cat# 179906
pHDM-Hgpm2 (Gag-Pol)	This study	Addgene Cat# 164441
pHDM-Tat1b (helper)	This study	Addgene Cat# 164442
pRC-CMV-Rev1b (helper)	This study	Addgene Cat# 164443
<b>Software and algorithms</b>		
GraphPad Prism 9.1.2	Graphpad Software	<a href="http://www.graphpad.com/scientific-software/prism/">www.graphpad.com/scientific-software/prism/</a> ; RRID:SCR_002798
Geneious Prime 2021	Geneious	<a href="http://www.geneious.com/">http://www.geneious.com/</a> ; RRID:SCR_010519
JMP Pro 15	SAS Institute	<a href="http://www.jmp.com/en_us/software/jmp.html">http://www.jmp.com/en_us/software/jmp.html</a> ; RRID:SCR_014242
FlowJo 10	FlowJo	<a href="https://www.flowjo.com/">https://www.flowjo.com/</a> ; RRID:SCR_008520
Fluent Control	Tecan	<a href="https://lifesciences.tecan.com/fluent-laboratory-automation-workstation?p=tab-4">https://lifesciences.tecan.com/fluent-laboratory-automation-workstation?p=tab-4</a>
R v4.0.2	Open source software	<a href="https://cran.r-project.org/bin/windows/base/">https://cran.r-project.org/bin/windows/base/</a> ; RRID:SCR_001905
CellCapTure	Stratedigm	<a href="https://stratedigm.com/cellcapture/">https://stratedigm.com/cellcapture/</a>
QuantStudio 12K Software v1.3	Applied Biosciences	<a href="https://www.thermofisher.com/us/en/home/global/forms/quantstudio-12k-flex-software-download.html">https://www.thermofisher.com/us/en/home/global/forms/quantstudio-12k-flex-software-download.html</a>

**RESOURCE AVAILABILITY**

**Lead contact**

Further information and requests for resources and reagents should be directed to, and will be fulfilled by, Alejandro Balazs ([abalazs@mgh.harvard.edu](mailto:abalazs@mgh.harvard.edu)).

**Materials availability**

All unique plasmids generated in this study are available through Addgene. Antibodies and other reagents are available from their respective sources.

**Data and code availability**

This study did not generate sequence data or code. Data generated in the current study (including neutralization and antibody measurements) have not been deposited in a public repository but are available from the lead contact upon request.

## EXPERIMENTAL MODEL AND SUBJECT DETAILS

### Human subjects

Use of human samples was approved by Partners Institutional Review Board (protocol 2020P002274). Serum samples from 239 COVID-19 vaccinees were collected. For each individual, basic demographic information including age and sex as well as any relevant COVID-19 history was obtained.

### Cell lines

HEK293T cells (ATCC) were cultured in DMEM (Corning) containing 10% fetal bovine serum (VWR), and penicillin/streptomycin (Corning) at 37°C/5% CO<sub>2</sub>. 293T-ACE2 cells were a gift from Michael Farzan (Scripps Florida) and Nir Hacohen (Broad Institute) and were cultured under the same conditions. Confirmation of ACE2 expression in 293T-ACE2 cells was done via flow cytometry.

## METHOD DETAILS

### Construction of variant spike expression plasmids

To create Delta and Omicron variant spike expression plasmids, we performed multiple PCR fragment amplifications utilizing oligonucleotides containing each desired mutation (Azenta) and utilized overlapping fragment assembly to generate the full complement of mutations for each strain. Importantly, we generated these mutations in the context of our previously described codon-optimized SARS-CoV-2 spike expression plasmid harboring a deletion of the C-terminal 18 amino acids that we previously demonstrated to result in higher pseudovirus titers. Assembled fragments were inserted into NotI/XbaI digested pTwist-CMV-BetaGlobin-WPRE-Neo vector utilizing the In-Fusion HD Cloning Kit (Takara). All resulting plasmid DNA utilized in the study was verified by whole-plasmid deep sequencing (Illumina) to confirm the presence of only the intended mutations.

### SARS-CoV-2 pseudovirus neutralization assay

To compare the neutralizing activity of vaccinee sera against coronaviruses, we produced lentiviral particles pseudotyped with different spike proteins as previously described (Garcia-Beltran et al., 2021b). Briefly, pseudoviruses were produced in 293T cells by PEI transfection of a lentiviral backbone encoding CMV-Luciferase-IRES-ZsGreen as well as lentiviral helper plasmids and each spike variant expression plasmid. Following collection and filtering, production was quantified by titrating via flow cytometry on 293T-ACE2 cells. Neutralization assays and readouts were performed on a Fluent Automated Workstation (Tecan) liquid handler using 384-well plates (Grenier). Three-fold serial dilutions ranging from 1:12 to 1:8,748 were performed for each serum sample before adding 50–250 infectious units of pseudovirus for 1 h. Subsequently, 293T-ACE2 cells containing polybrene were added to each well and incubated at 37°C/5% CO<sub>2</sub> for 48 h. Following transduction, cells were lysed using a luciferin-containing buffer (Siebring-van Olst et al., 2013) and shaken for 5 min prior to quantitation of luciferase expression within 1 h of buffer addition using a Spectramax L luminometer (Molecular Devices). Percent neutralization was determined by subtracting background luminescence measured in cell control wells (cells only) from sample wells, and dividing by virus control wells (virus and cells only). Data were analyzed using Graphpad Prism, and NT50 values were calculated by taking the inverse of the 50% inhibitory concentration value for all samples with a neutralization value of 80% or higher at the highest concentration of serum. NT50 values were converted to WHO International Units using the Human SARS-CoV-2 Serology Standard (Lot #COVID-NS01097) obtained from NCI-Frederick National Laboratory for cancer research, which was calibrated to the WHO SARS-CoV-2 Serology International Standard (20/136).

### Titering and flow cytometry

To determine the infectious units and infection rates of SARS-CoV-2 pseudovirus, we plated 400,000 293T-ACE2 cells per well of a 12-well plate. SARS-CoV-2 pseudovirus supernatants that were produced in parallel under identical conditions (as described above) were added 24 h later in three ten-fold serial dilutions. Cells were incubated for 48 h at 37°C/5% CO<sub>2</sub> to allow for expression of ZsGreen reporter gene and harvested with Trypsin-EDTA (Corning). Cells were resuspended in PBS supplemented with 2% FBS (PBS+), and analyzed on a Stratadigm S1300Exi Flow Cytometer to determine the percentage of ZsGreen<sup>+</sup> cells. Infectious units were calculated by determining the percentage of infected cells in wells exhibiting linear decreases in transduction and multiplying by the average number of cells per well determined at the initiation of the assay. At low MOI (infection rates <10%), each transduced ZsGreen<sup>+</sup> cell was assumed to represent a single infectious unit. Surface expression of ACE2 on 293T and 293T-ACE2 cells was confirmed by staining with 1 μL of anti-human ACE2, Alexa Fluor 647-conjugated antibody (R&D). Cells were incubated for 10 min at 25°C prior to running on a Stratadigm S1300Exi Flow Cytometer.

### Quantitation of pseudovirus by RT-qPCR

To determine the genome copy concentration of pseudotyped lentiviral vectors, lentiviral RNA was extracted from pseudovirus supernatant using the QIAamp viral RNA mini kit (Qiagen). Each sample was serially diluted, and each dilution was treated with 1.2 U of Turbo DNase (Invitrogen) at 37°C for 30 min followed by heat inactivation at 75°C for 15 min. 10 μL of the treated RNA was used in a 20 μL qRT-PCR reaction with the qScript XLT one-step RT-qPCR Tough Mix, low ROX mix (Quanta Biosciences), a TaqMan probe containing locked nucleic acids (56-FAM/AGC+C/i5NiIInd/GG+GA/ZEN/GCTCTCTGGC/3IABkFQ/) (IDT), and primers designed for

targeting the LTR gene of NL4-3 HIV genome, from which the lentiviral vector was derived (5'-GGTCTCTCTIGITAGACCAG and 3'-TTTATTGAGGCTTAAGCAGTGGG). Each dilution was run in duplicate on a QuantStudio 12K Flex (Applied Biosystems). The following cycling conditions were used: 50°C for 10 min, 95°C for 3 min followed by 50 cycles of 95°C for 15 s and 60°C for 1 min. Virus titer was determined by comparison with a standard curve generated using a plasmid standard generated from serial dilution of CMV-Luciferase-IRES-ZsGreen lentiviral backbone. DNase and No DNase controls were also included at  $2.5 \times 10^8$  GC/mL of the same plasmid. The range of the assay was from  $2.5 \times 10^7$  GC/mL to  $1.5 \times 10^3$  GC/mL. Upon analysis, the average of the three most concentrated dilutions within range of the standard were used to calculate genome copies/mL.

#### **Anti-SARS-CoV-2 spike and nucleocapsid antibody assays**

The EUA-approved electrochemiluminescence-based Roche Elecsys® anti-SARS-CoV-2 spike antigen (semi-quantitative) and nucleocapsid antigen (qualitative) immunoassays were used to detect total antibodies (IgG, IgM, and/or IgA antibodies) to SARS-CoV-2 spike and nucleocapsid in vaccinee sera. The upper limit of detection for the spike semi-quantitative assay was 25,000 U/mL. The assay was run on a Roche Cobas 8000 e801 Immunoassay Analyzer in the Massachusetts General Hospital Core Laboratory.

#### **QUANTIFICATION AND STATISTICAL ANALYSIS**

Data and statistical analyses were performed using GraphPad Prism 9.2.0, JMP Pro 16.1.0 (SAS Institute), and R v4.0.2. Flow cytometry data were analyzed using FlowJo 10.7.1. Student's *t* tests were performed in the indicated figures; all *p* values are unadjusted. Statistical significance was defined as  $p < 0.05$ . Error bars throughout all figures represent one standard deviation or standard error of the mean where indicated.